

## Procedure for saliva sample collection:



### Collection precautions:

Donor should not eat, drink, smoke or chew gum for 30 minutes before collecting saliva sample.

Ensure the sponge tip does NOT come into contact with any surface prior to collection.

**Intended use:** ORACollect®-Dx is intended for use in the non-invasive collection of saliva samples. Human DNA from the saliva sample is isolated, stabilized, and suitable for use in FDA cleared molecular diagnostic applications. Saliva samples collected using ORACollect-Dx are stabilized and can be transported and/or stored long-term at ambient conditions.

**Contents:** Kit contains stabilizing liquid.


### Warning and precautions:

- Rx only. For professional use.
- Read all instructions carefully prior to sample collection; deviation may result in inadequate sample and impact DNA yield.
- Choking hazard. Caution should be used when inserting sponge into the mouth.
- Donors with xerostomia (dry mouth) may not collect adequate sample using these instructions resulting in lower DNA yield and an invalid sample for use.
- Wash with water if stabilizing liquid comes in contact with eyes or skin. Do NOT ingest. See MSDS at [www.dnagenetek.com](http://www.dnagenetek.com)


**Storage:** 15-25°C

### Label legend:


	Catalog number
	In vitro diagnostic medical device
	CE Marking
	Manufacturer
	Authorized Representative
	Storage instructions
	Collect sample by (Use by)
	Lot number

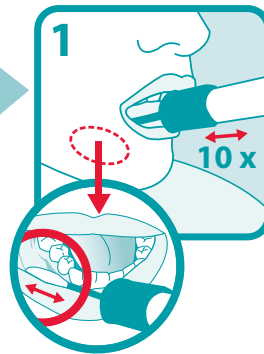


Open package and remove collector without touching sponge tip. Place sponge as far back in the mouth as comfortable and rub along the lower gums (see close up image) in a back and forth motion. Gently rub the gums 10 times. If possible, avoid rubbing the teeth.

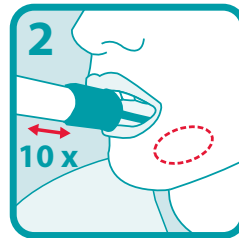


Gently repeat rubbing motion on the opposite side of the mouth along the lower gums for an additional 10 times.

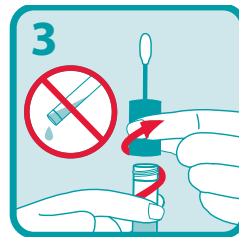

Hold the tube upright to prevent the stabilizing liquid inside the tube from spilling. Unscrew the blue cap from the collection tube without touching the sponge.

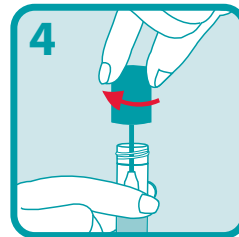

Turn the cap upside down, insert the sponge into the tube and close cap tightly.

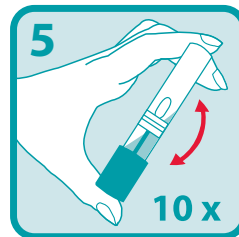

Invert the capped tube and shake vigorously 10 times.



For In Vitro Diagnostic Use

**DNAgenetek**



Made in Canada  
DNA Genetek Inc.  
3000 - 500 Palladium Drive  
Ottawa, ON, Canada K2V 1C2

*Superior samples*  
*Proven performance*

Toll-free (North America): 1.866.813.6354  
Tel: +1.613.723.5757 • Fax: +1.613.723.5057  
info@dnagenetek.com  
www.dnagenetek.com



Emergo Europe, Prinsessegracht 20, 2514 AP The Hague, The Netherlands  
Australian Sponsor: Emergo Australia, Level 20, Tower II, Darling Park, 201 Sussex Street, Sydney, NSW 2000 Australia

\*ORACollect is a registered trademark of DNA Genetek Inc.  
Some DNA Genetek products may not be available in all geographic regions, contact your sales representative for details.  
All DNA Genetek protocols, white papers and application notes, are available in the support section of our website at [www.dnagenetek.com](http://www.dnagenetek.com).

Patent ([www.dnagenetek.com/legalnotices](http://www.dnagenetek.com/legalnotices))  
© 2017 DNA Genetek Inc., a subsidiary of OraSure Technologies, Inc., all rights reserved.  
PD-PR-00506 Issue 4/2017-05