
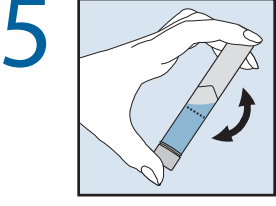
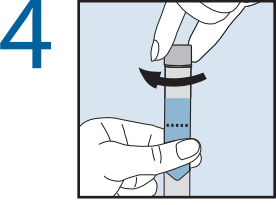
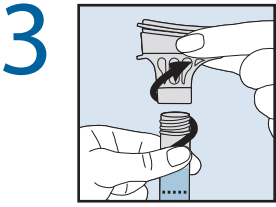
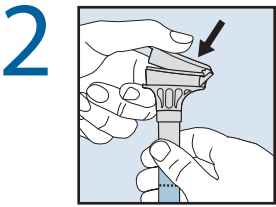
For Research Use Only
Not for use in diagnostic procedures

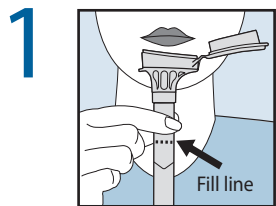
DNAGENOTEK™

■ DNA Genotek Inc.
3000 - 500 Palladium Drive
Ottawa, ON, Canada K2V 1C2

Tel.: +1.613.723.5757
info@dnagenotek.com
www.dnagenotek.com

Patent (www.dnagenotek.com/
legal-notice)
PD-PR-617 v2

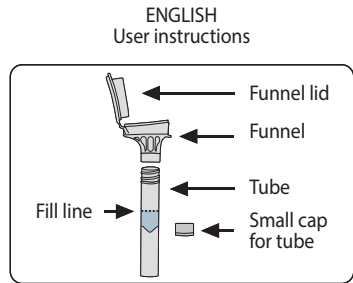**Read all instructions
prior to collection.**

Intended use: For the collection of human DNA from saliva samples.

Contents: 1 kit containing 2 mL stabilizing liquid.

Warnings and precautions:

Wash with water if stabilizing liquid comes in contact with eyes or skin. Do NOT ingest. See SDS at www.dnagenotek.com.

Small cap may pose a choking hazard.

Storage: 15°C / 30°C (59°F / 86°F)

Summary and explanation of the kit:

Oragene™DISCOVER is a self-collection kit that provides the materials and instructions for collecting and stabilizing saliva specimens.

Collection precautions:

Do NOT eat, drink, smoke, or chew gum for 30 minutes before giving your saliva sample.

Do NOT remove the plastic film from the funnel lid.


Procedure:

Most people take between 2 and 5 minutes to deliver a saliva sample following steps 1 to 5.

- 1 Spit into funnel until the amount of liquid saliva (not bubbles) reaches the fill line shown in picture #1.
- 2 Hold the tube upright with one hand. Close the funnel lid with the other hand (as shown) by firmly pushing the lid until you hear a loud click. The liquid in the lid will be released into the tube to mix with the saliva. Make sure that the lid is closed tightly.
- 3 Hold the tube upright. Unscrew the funnel from the tube.
- 4 Use the small cap to close the tube tightly.
- 5 Shake the capped tube for 5 seconds. Discard the funnel.

Label legend:

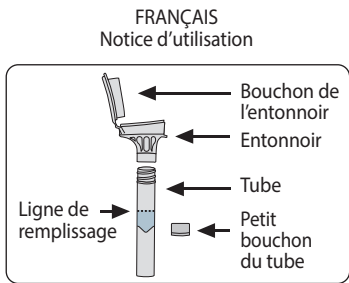- Flammable liquid
- Consult package insert
- Collect saliva by (Use by)
- Catalog number
- Caution, consult instructions for use
- 15°C / 30°C Storage instructions
- 59°F / 86°F
- Manufacturer
- Do not re-use
- Lot number

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all rights reserved.