

### Kakaretso le tlhaloso ya disebediswa:

OMNIGENE-GUT e tswelletsa dihlahiswa le ditaello bakeng sa ho bokella le ho lekalekanya microbial DNA bakeng sa disampole tsa mantle.

### Ditemoso le ditlhokomediso:

- HA E NOWE, E SEBEDISWA MMELENG.
- O SE KE wa tlosa tube e tshehla ya ka hodimo ho tube.
- O SE KE wa tsholla metsi a ho lekalekanya a tube.
- Hlatswa ka metsi haeba metsi a sesebediswa a ka kgoma mahlo kapa letlalo. O SE KE wa nwa.
- Haeba o bokella mohlala wa mantle a leng jwalo ka metsi, sheba ditlhaloso tse rometsweng tsa ka thoko.
- Dihlahiswa tse dinyenyane di ka etsa hore o bethwe ke motlakase.

**Bobolokelo:** 15°C ho ya ho 25°C

Romela ho ya ka melao ya sebaka ya ho romela ditlhalobo tsa sebaoloji ka dinamelwa. Sheba MSDS ho [www.dnagenotek.com](http://www.dnagenotek.com)

### Histori ya Label:

- Sampole e bokelletswa ke (Sebediswa ke)
- Nomoro ya khathaloko
- Moetsi
- Ditaello tsa ho boloka
- Temoso, sheba ditaello tsa tshebediso
- Nomoro ya Lot

## DITAELO TSA MOSEBEDISI

### Bala ditaello kaofela pele o bokella

#### Tsamaiso:



#### 1 DITOKISETSO TSE BOHLOKWA:

- O lokela ho rota pele o ka qala ho bokella.
- Bokella mohlala wa mantle a se nang moroto kapa metsi a ntlwaneng ya ho ithusa.
- Pampiri ya ntlwana kapa thishu di ka sebediswa.




- #### 2
- Ha o ntse o tshwere bohodimo ba tube e tshehla, kurufulla FEELA khepe e phepholo mme o e bee bakeng sa ho sebedisa nakong e tlang.



#### TSE BOHLOKWA:

O SE tlose bohodimo bo tshehla ba tube  
O SE tsholle metsi a tekatekanyo a tube.

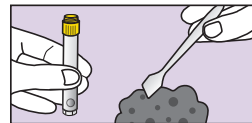

- #### 3
- Sebedisa spatula ho bokella sampole e nyenyane ya mantle.



Boholo ba sampole ya mantle.



- #### 4
- Fetsetsa sampole ya mantle ka hodimo ha tube e tshehla. Pheta sena ho fihlella sampole e tlatsa bohodimo ba tube e tshehla.



**SE BOHLOKWA:** O SE ke wa leka ho kgarametsa sampole ka ho tube.



- #### 5
- Hlophisa ka ho rapameng ho bohodimo ba tube hore sampole e lekalekane mme o tlose e feteletseng.

Phumula bokantle ba bohodimo ba tube ka pampiri ya ntlwana kapa thishu ha ho hlokahala.




- #### 6
- Nka khepe e phepholo e nang le qetello e tiileng e shebileng tlase mme o e kurufele ho bohodimo ba tube e tshehla ho fihlella e kwetse hantle.



- #### 7
- Tsukutla tube e kwetsweng ka thata ka mokgwa wa pele le morao ka metsotswana e 30.

Metsi e 30.

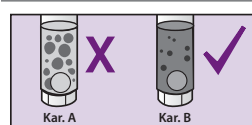

- #### 8
- Sampole ya mantle e tla kopana le metsi a tekatekano ho tube; ha se dintho kaofela tse tla qhibilisa.

**TSE BOHLOKWA:** Tswelapele ho tsukutla haeba dikarolo tse **kgolo** di ka ba di sa le teng ho ya ka karolo A.



- #### 9
- Bea spatula ho sephuthelwana sa yona sa ho qala kapa se dikolose ka pampiri ya ntlwana mme o lahle ho ditshila.

**SE BOHLOKWA:** Romela sampole hore e shebiwe o latela ditaello tsa ho romela tse hlalishitsweng ho kit ya hao.